

Optical Temperature Sensors

USER MANUAL

OPERATING INSTRUCTIONS



PyroScience optical temperature sensors are manufactured by

PyroScience GmbH
Hubertusstr. 35
52064 Aachen
Germany

Phone	+49 (0)241 5183 2210
Fax	+49 (0)241 5183 2299
Email	info@pyroscience.com
Internet	www.pyroscience.com

Registered: Aachen HRB 17329, Germany

TABLE OF CONTENT

1	INTRODUCTION	4
2	QUICK START	5
3	SENSOR SETTINGS.....	6
4	SENSOR CALIBRATION	7
4.1	CALIBRATION PROCEDURE	8
5	SENSOR APPLICATION.....	9
5.1.1	<i>Fiber-based sensors</i>	<i>9</i>
5.1.2	<i>Contactless and combined sensors</i>	<i>10</i>
6	CLEANING AND STORAGE.....	11
6.1	STERILIZATION	11
6.2	CLEANING AND STORAGE	11
7	RELATED DOCUMENTS	12
8	APPENDIX.....	13
8.1	TEMPERATURE MEASURING PRINCIPLE	13
8.2	EXPLANATION OF THE SENSOR CODE	13
8.3	AVAILABLE SENSORS AND READ-OUT DEVICES	15
8.4	PT100 TEMPERATURE SENSOR CALIBRATION	16
9	WARNINGS AND SAFETY GUIDELINES	17

1 Introduction

PyroScience offers fiber-based and contactless optical temperature sensors, as well as a number of combined sensors like flow-through cells respiration vials and sensor vials for simultaneous measurements of temperature, oxygen and pH.

<https://www.pyroscience.com/en/products/all-sensors/attributes/Temp>

These sensors can be read-out with the following fiber-optic meters from *PyroScience*

- the multi-channel PC-operated *FireStingO2* (with *Pyro Oxygen Logger* software) for optical O₂ and temperature sensors and
- the multi-analyte & multi-channel PC-operated *FireSting pro* (with *Pyro Workbench*) for optical pH, O₂ and temperature sensors.

All software versions are available for free download from the *PyroScience* website and must be installed on a Windows PC/laptop before connecting the respective meter for the first time. For details on the read-out devices and their software, please see their respective manuals and handling guidelines.

This manual is intended to provide all necessary information on standard application of optical temperature sensors from *PyroScience*.

For more information concerning advanced applications, please contact us at info@pyroscience.com.

Your *PyroScience* Team

2 Quick Start

Step 1: Please download the correct software in the downloads tab of your purchased device on www.pyroscience.com, unzip and start the installer, and follow the instructions.

Step 2: After successful installation, connect the *PyroScience* read-out device to the Windows PC/laptop with the micro-USB cable.

Step 3: Carefully remove the protective caps from the sensor tip, from the fiber plug and from the optical connector(s) at the *FireSting* read-out device.

Step 4: Connect the *PyroScience* temperature sensor(s) to the **optical channel connector(s)** of the device.

Step 5: For **calibration of optical temperature sensors**, connect an appropriate Pt100 temperature sensor to the temperature port. Alternatively, a fixed temperature can be entered (needs to be measured and kept constant).

Step 6: Start the *PyroScience* software by clicking on the short-cut on your desktop.

Step 7: Enter all required **Sensor Settings** for each connected sensor, including the **sensor code**, **sample interval**, **measuring mode** and, for contactless sensors, also the **fiber length (m)**.


Step 8: Perform a 1-point **Calibration** of the optical temperature sensor. As reference temperature, select either the connected Pt100 temperature sensor (needs to be calibrated first, see chapter 8.4) or enter a fixed temperature at constant temperature conditions (e.g. in a water bath).

Step 9: Start measurements and activate **Data Logging**.

3 Sensor Settings

Each optical temperature sensor comes with an individual **Sensor Code**, containing important information for optimal sensor settings and for calibration. The first letter of the sensor code defines the sensor type. Therefore, it is important to enter the Sensor Code of the connected sensor into the **Sensor Settings** of the respective software. The number of the channel tab must correspond with the channel number at the *FireSting* read-out device.

IMPORTANT: Enter the correct **Sensor Code** for sensors connected to a channel at a *FireSting* read-out device. The sensor code can be found on the label attached to the cable (fiber-based sensors) or on the bag of contactless sensors (see example below).

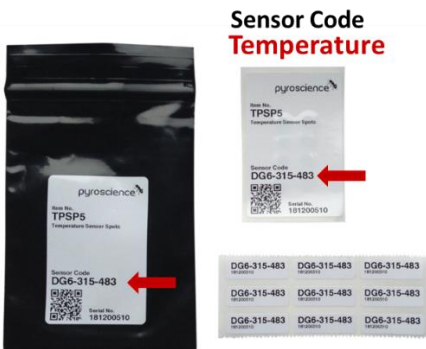

Next the **Measuring Mode** needs to be entered (refer to the *Pyro Workbench* or *FireStingO₂* manual for further information). This can be adjusted gradually between low drift and low noise of the sensor signal by moving the arrow with the mouse along the scale. An intermediate mode is default.

The **Unit** of the temperature readings is **°C**.

4 Sensor Calibration

Ensure that the correct **Sensor Code** has been entered in the settings (refer to chapter 3).

The 1-point calibration of optical temperature sensors can be performed against two different reference temperatures:

- the temperature measured by the external **Pt100 Temperature Sensor** (item *TDIP15* or *TSUB21*)
- a **Fixed Temperature** that needs to be determined and entered in the software (must be kept constant)

For precise optical temperature sensor calibration using the External Temperature Sensor, it needs to be determined manually if the external (Pt100) temperature sensor has an offset. In case of an offset, the Pt100 temperature sensor needs to be calibrated first (see Appendix 8.4) before calibrating the optical sensor.

If a Manual Temperature was selected, the temperature in the sample/calibration standard must be measured, adjusted and kept constant (needs to be controlled)! Ensure constant and defined conditions!

- Gas measurements (only in combination with oxygen sensors): the sensor should be calibrated in ambient air at defined temperature
- Measurements in aqueous/water samples: the sensor has to be calibrated in water at defined temperature

NOTE: It is strongly recommended to perform the calibration at conditions close to the environmental temperature during measurements. Ensure constant temperature during calibration!

4.1 Calibration procedure

NOTE: The calibration of the optical temperature sensor (recommended every day) is only possible if the **correct Sensor Code** has been entered in the **Settings**.

Calibration should be performed following the instructions of the software (*Pyro Workbench*) or *FireStingO2* manual (*Oxygen Logger Software*). It is generally recommended to perform a 1-point calibration in gas (*water*) for gas (*water*) measurements.

Important: The device and sensors must be placed **for >30 min. under constant temperature conditions** before the calibration is performed.

For a calibration with the external temperature sensor (item no. *TDIP15* or *TSUB21*), ensure that the external temperature sensor is placed close to and exactly at the same temperature conditions as the optical temperature sensor.

For calibration against a fixed/manual temperature (°C), use a calibrated incubator (for measurements in the gas phase) or water bath of constant and known temperature (for measurements in water/aqueous samples).

When the sensor is placed into the calibration medium (air or water), wait until the sensor reading is stable by observing the graph and the numerical display of the temperature sensor reading. If all readings have reached their steady-state, calibrate the optical temperature sensor accordingly.

A completed calibration is indicated by the green indicator.

5 Sensor Application

PyroScience optical temperature sensors can be applied in gas phases and aqueous solutions, and are typically used for two different purposes:

- for precise and fast temperature measurements in water
- for true temperature compensation of *PyroScience* optical oxygen and pH sensor readings

Specific application and calibration instructions for different sensors are listed in the following table.

5.1.1 Fiber-based sensors

Sensor item	Sensor-specific application instructions
<i>TPR430</i>	Application: water & gas & semi-solid samples Calibration: 1-point calibration Note: Handle with care to avoid injuries! Extend sensor tip for calibration and measurements.
<i>TPF1100</i>	Application: water & gas & semi-solid samples (only in combination with item <i>OXF1100</i>) Calibration: 1-point calibration Note: Handle with care! Unprotected fragile sensor tip.
<i>TPB430</i>	Application: water & gas & custom application Calibration: 1-point calibration Note: Handle with care, especially during custom integration! Unprotected fragile sensor tip. Avoid breakage!

5.1.2 Contactless and combined sensors

Sensor item	Sensor-specific application instructions
<i>TPSP₅</i>	<p>Application: water & gas (gas only with item <i>OXSP₅</i>)</p> <p>Calibration: 1-point calibration</p> <p>Features: optical isolation</p> <p>Sterilization: 70% ethanol, 70% isopropanol, can be autoclaved (details on request)</p> <p>Note: Mind air bubbles! Glue carefully with silicone glue and let dry for 24h.</p>
<i>TOVIAL₂₀</i> <i>TOVIAL₄</i>	<p>Application: water & gas</p> <p>Calibration: 1-point for temperature sensor, 1- or 2-point for oxygen sensor (see oxygen sensor manual)</p> <p>Features: optical isolation, O₂ & Temp</p> <p>Sterilization: 70% ethanol, 70% isopropanol</p> <p>Note: Remove air bubbles! Determine specific volume before measurements. Ensure stable temperature conditions.</p>
<i>TOFTC₂</i>	<p>Application: water & gas</p> <p>Calibration: 1-point for temperature sensor, 1- or 2-point for oxygen sensor (see oxygen sensor manual)</p> <p>Features: optical isolation, O₂ & Temp</p> <p>Sterilization: 70% ethanol, 70% isopropanol</p> <p>Note: Liquid flow rate of 20-500 mL/min recommended. Remove air bubbles! Clean regularly.</p>
<i>PHTOVIAL₂₀</i>	<p>Application: water</p> <p>Calibration: 1-point for temperature sensor, 1- or 2-point in specified pH buffers/oxygen calibration standards (see pH and oxygen sensor manual)</p> <p>Features: optical isolation, pH & O₂ & Temp</p> <p>Note: Remove air bubbles! Determine specific volume before measurements. Ensure stable temperature conditions.</p>

6 Cleaning and Storage

6.1 Sterilization

Refer to the specifications on the respective *PyroScience* website.

The temperature sensor spots (item no. *TPSP5*) can be autoclaved. More details on request. The temperature sensor spots and the combined sensors can be treated with 70% ethanol and 70% isopropanol.

IMPORTANT: do **not** use bleach, acetone or any solvent/agent not approved by *PyroScience*!

6.2 Cleaning and Storage

After the measurements, the sensor tip of the needle-type and bare fiber sensors should be rinsed carefully with demineralized water. The sensors can further be cleaned with 3% H₂O₂ and soap solution, followed by thorough rinsing with demineralized water.

After cleaning, let dry and put on the protective cap/tubing for storage in a dry, dark and secure place at room temperature. Put the black caps on the plug of the fiber to prevent that light is entering the fiber possibly causing photo-bleaching of the indicator.

For retractable sensors and application in samples with dissolved salts (e.g. seawater), the sensor must be cleaned thoroughly with demineralized water to prevent salt crystallization in the needle which can cause breaking of the sensor tip. **After drying**, retract the sensor tip into the needle and put on the protective cap onto the needle to protect the sensor tip and to avoid injuries.

Store the sensor in a dry, dark and secure place at room temperature.

7 Related documents

Related documents for more detailed instructions on fiber-optic read-out devices, software and optical sensors are available:

- manual for logger software "*Pyro Workbench*" (Windows)
- manual for multi-analyte meter *FireSting pro*
- manual for oxygen meter *FireStingO₂* (with *Oxygen Logger* software)
- manual for optical pH sensors
- manual for optical oxygen sensors

8 Appendix

8.1 Temperature Measuring Principle


Optical temperature sensors rely on thermal quenching of photoluminescence. In this physical process, luminescent properties of the material are reversibly affected by temperature. The sensor response times are only limited by the dimensions of the sensor and are ultra-fast (<0.1 s in water), like in case of fiber-optic minisensors.

The utilized temperature-sensitive material is specially optimized to be fully compatible with the innovative *REDFLASH* technology from *PyroScience*, being excitable with orange red light (610-630 nm) and showing emission in the near infrared part of the spectrum. Importantly, the optical temperature sensors are characterized by extreme photo-stability, which enables virtually unlimited number of measurements, and does not show interferences with gases (such as oxygen) and ionic species.

8.2 Explanation of the Sensor Code

The temperature sensors are delivered with an attached sensor code which must be entered into the Settings (refer to chapter 3). The following figure gives a short explanation of the information given in the sensor code.

Example Code: **DG6-310-479**

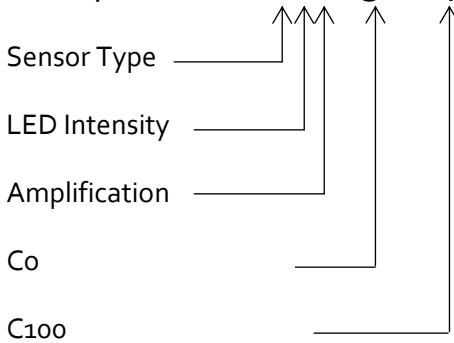

Sensor Type

- C Optical Temperature Minisensor (*Thermoblue*)
- D Optical Temperature Sensor Spot/FTC (*Thermogreen*)

LED Intensity

- | | | | |
|---|-----|---|------|
| A | 10% | E | 40% |
| B | 15% | F | 60% |
| C | 20% | G | 80% |
| D | 30% | H | 100% |

Amplification

- | | |
|---|------|
| 4 | 40X |
| 5 | 80X |
| 6 | 200X |
| 7 | 400X |

Optical Temperature Sensors

Co and **C100**: values needed for calibration.

8.3 Available sensors and read-out devices

Sensor Type	Compatible read-out devices				
	<i>FSO2-X</i>	<i>FSGO2</i>	<i>FSO2-SUBPORT*</i>	<i>Piccolo2</i>	<i>FS pro</i>
<i>TPR430</i>	✓	✗	✓	✗	✓
<i>TPF1100</i>	✓	✗	✓	✗	✓
<i>TPB430</i>	✓	✗	✓	✗	✓
<i>TPSP5</i>	✓	✗	✓	✗	✓
<i>TOFTC2</i>	✓	✗	✗	✗	✓
<i>TOVIAL4</i> <i>TOVIAL20</i>	✓	✗	✗	✗	✓
<i>PHTOVIAL20</i>	✓	✗	✗	✗	✓

*only in combination with UNDERWATER CONNECTOR (option - SUB)

8.4 Pt100 Temperature Sensor Calibration

For precise absolute temperature readings and precise calibration of optical temperature sensors, an optional **1-point calibration of the external temperature sensor Pt100** is recommended.

For this, check the reading of the external temperature Pt100 probe periodically in stirred water / water bath / incubator of known temperature at steady state. It is also possible to prepare a water-ice-mixture giving 0°C, where at least 50 mm of the Pt100 temperature probe tip is submerged.

IMPORTANT: After calibration of the Pt100 temperature sensor, a new optical sensor calibration must be performed.

9 Warnings and Safety Guidelines

Before using *PyroScience* temperature sensors, carefully read the instructions and user manuals for the respective *PyroScience* read-out device. Please find the correct device manual in the downloads tab of your purchased device on www.pyroscience.com

Prevent mechanical stress (e.g. scratching) to the sensing surface at the tip of the temperature sensor! Avoid strong bending of the fiber-optic cables. They might break!

Ensure that the complete sensing surface at the tip is always covered by the sample and is free of air bubbles, and that liquid samples are stirred.

Calibration and application of the temperature sensors is on the user's authority, as well as data acquisition, treatment and publication!

PyroScience temperature sensors and read-out devices are not intended for medical or military purposes or any other safety-critical applications. They must not be used for applications in humans; not for in vivo examination on humans, not for human-diagnostic or therapeutic purposes. The sensors must not be brought in direct contact with foods intended for consumption by humans.

The sensors must be used in the laboratory by qualified personnel only, following the user instructions and the safety guidelines of the manual, as well as the appropriate laws and guidelines for safety in the laboratory!

Keep the *PyroScience* temperature sensors and read-out devices out of reach of children! Store the temperature sensors in a secure, dry and dark place at room temperature.